WIDEX **SMARTRIC**™

QUICK REFERENCE GUIDE

Welcome to **new angles of natural hearing** with the new Widex SmartRIC™ – a hearing aid with a smart design that further supports the hearing aid user audiologically.



This quick reference guide provides you with an overview of the most important information on the Widex SmartRIC hearing aids.

RECEIVERS AND EARWARE

The Widex SmartRIC hearing aid is compatible with the **V.2 M- and P-Receivers**.


It can be used with the following instant ear-tips:



Open



Open Round

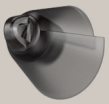

Tulip




Round



Sleeve Vented



Double

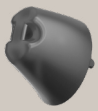

Sleeve Power

What about custom ear-tips?
Every ear (canal) is unique and that impacts the size of a custom ear-tip. It's important that you check the sizing versus space available in the charger.

ACCESSORIES

The Widex SmartRIC is also compatible with our comprehensive range of accessories.

Widex Sound Assist
Bluetooth streaming, remote control, table mic, telecoil functionality



Widex SoundConnect
Audio streaming from PCs and laptops



Widex TV Play
Audio streaming from televisions



Widex RC-DEX
Remote control – volume, program changes



APPS

Several Widex apps are also compatible with the Widex SmartRIC. These include:

WIDEX MOMENT™
Control volume and change programs; tailor sound with MySound

WIDEX REMOTE CARE™
Perform remote hearing aid fittings and adjustments

WIDEX TONELINK™
Basic remote control functions; for the easygoing user

LED LIGHT INDICATOR

The Widex SmartRIC comes with a nifty LED light indicator to show the pairing, battery and power status.

Light indicators	Status
Five slow green flashes	Hearing aid is turning on.
Six green flashes	Hearing aid has been paired to a Bluetooth device.
Continuous green flashes*	Hearing aid is turned on, battery > 20%.
Continuous red flashes*	Hearing aid is turned on, battery < 20%.

* Additional option that can be enabled.

FEATURES

The Widex SmartRIC is available in the 440, 330, 220 and 110 price points.

KEY FEATURES	440	330	220	110
Performance	●●●●●●	●●●●●	●●●●	●●
Processing & fine-tuning channels	15	12	10	6
Widex PureSound™ with ZeroDelay™ technology	●	●	●	●
TruAcoustics™ with Personal Gain Integrator	●	●	●	●

CONNECTIVITY	440	330	220	110
2.4 GHz control (Android and iOS)	●	●	●	●
2.4 GHz audio streaming (iOS and ASHA)	●	●	●	●
WidexLink to DEX assistive listening devices*	●	●	●	●

APPs FOR iOS & ANDROID	440	330	220	110
MOMENT App with SoundSense Learn	●	●	●	
MOMENT App	●	●	●	●
TONELINK App	●	●	●	●
COM-DEX App	●	●	●	●
Widex Remote Care App	●	●	●	●

FEATURES	440	330	220	110
Adaptation Manager	●	●	●	●
Fluid Sound Analyzer (Sound Classes)	11	7	4	3
Programs	5	4	3	3
Smartwind Manager	●			
High-frequency boost	●			
Speech Enhancer RT	RT/IE	IE		
HD Locator	●	●	●	
TruSound Softener	●	●	●	
SoundSense Adapt	●	●	●	
Preference control	●	●	●	●
LED light indicator	●	●	●	●
Soft-level noise reduction	●	●	●	●
Noise reduction	●	●	●	●
Zen/SoundRelax	●	●	●	●
Audibility Extender	●	●	●	●

*Incl. Widex Sound Assist, TV-PLAY, RC-DEX, TV-DEX, COM-DEX, FM+ DEX.